PCT

(30) Priority Data:

08/052,354

63134 (US).

WORLD INTELLECTUAL PROPERTY ORGANIZATION International Bureau



INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

THISTONAL WITHOUTS		
(51) International Patent Classification ⁵ :		(11) Interactional Publication Number: WO 94/25097
A61M 29/00	Al	(43) International Publication Date: 10 November 1994 (10.11.94)
(21) International Application Number: PCT/US	94/045	74 (81) Designated States: AU. CA. IP. European parent (AT. BE. CH, DE, DK, ES, FR, GB, GR, IE, IT, LU, MC, NL, PT.
(22) International Filips Date: 26 April 1994	(26.04.9	

US

(71) Applicant: MALLINCKRODT MEDICAL, INC. [US/US]; 675 McDonnell Boulevard, P.O. Box 5840, St. Louis, MO

26 April 1993 (26.04.93)


(72) Inventors: BARCLAY, Anita, L.; 1208 East Kyle Road, Clute, TX 77531 (US). LOBDELL, Robert; 5 Isabella Court, Angleton, TX 77515 (US).

(74) Agents: HEY, David, A. et al.; Mallinckrodt Medical, Inc., 675 McDonnell Boulevard, P.O. Box 5840, St. Louis, MO 63134 (US).

Published

With international search report. Before the expiration of the time limit for amending the claims and to be republished in the event of the receipt of amendments.

(54) Title: METHOD OF JOINING PORTIONS OF INFLATION CATHETER

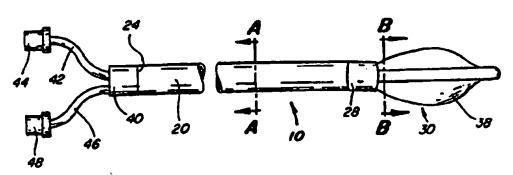

(57) Abstract

The present invention relates to a method of joining portions of an inflation/expansion catheter (10). In particular, the present invention relates to an improved method of joining coaxial lumen portions (35, 37) to multi-humen portions (50, 53) of an inflation/expansion catheter. More particularly, the present invention relates to a method of joining a multi-lumen shaft to a coaxial balloon shaft of a dilatation catheter.

FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AT	Austria	GB	United Kingdom	MR	Marritania
ΑÜ	Ameralia	GE	Georgia	MW	Mahwi
BB	Barbados	GN	Guinea	NE	Niger
BE	Belgium	GR	Groupe	NL	Netherlands
IJ	Burkina Paro	ETU	Hungary	NO	Norway
BG	Bulgaria	IR.	Ireland	NZ	New Zesland
Ŋ	Bento	<u> </u>	Italy	PL	Poland
BR	Brazil	JP	Japan	PT	Portugal
BY	Betarus	KE	Kenya	RO	Romanta
CA	Canada	KG	Kyrgystan	RU	Russian Federation
G₹	Central African Republic	KP	Democratic People's Republic	SD	Sodan
CG	Congo	_	of Korea	SE	Sweden
CH CH	Switzerland	KOR	Republic of Korea	SI	Stovenia
	Côte d'ivoire	KZ	Kazathstan	SK	Slovakia
CI CM		ũ	Liechennein	SN	Smept
	Camuroon	LK	Sri Lanks	TD	Chard
CN .	China	ᇤ	Lexembourn	TG	Togo
ß	Czechodovskia	LV	Latvia	IJ	Tajikistan
CZ	Czech Republic			π	Trinidad and Tobago
DE	Germany	MC	Monaco	ÜA	Utraine
DK	Deemark	MD	Republic of Moldova		
ES	Spain	MG	Madegracer	C/B	United States of America
n	Pintand	MIL	Mail	UZ	Uzbekistan
т	France	MIN	Mongolia	VN	Vist Nam
GA	Gebon				

WO 94/25097 PCT/US94/04574

METHOD OF JOINING PORTIONS OF INFLATION CATHETER

Background

The present invention relates to a method of joining portions of an inflation/expansion catheter. In particular, the present invention relates to an improved method of joining coaxial lumen portions to multi-lumen portions of an inflation/expansion catheter. Further, the present invention relates to a dilatation catheter formed by such method.

transluminal angioplasty (PTA) Percutaneous procedures typically are performed by introducing a dilatation balloon catheter into the cardiovascular system of the patient through the brachial or femoral arteries via a sheath introducer. In particular, a guidewire and the dilatation balloon catheter are introduced through the sheath introducer, the guidewire being disposed within an inner lumen of the balloon The guidewire and balloon are advanced until catheter. the dilatation balloon is properly located within the area of lesion to be treated. Once positioned, the inflated is pressure dilatation balloon predetermined size using a radiopaque liquid, such as contrast medium, through a second lumen of the balloon radially compress order to in catheter. atherosclerotic plaque of the lesion against the inside of the artery wall and thereby dilate the lumen of the The dilatation balloon is then deflated. Inflation and deflation of the balloon may be repeated several times within the lesion until the desired results The balloon catheter is then removed so are achieved. that blood flow may be resumed through the dilated artery.

It is desirable to provide a method of joining coaxial lumen and multi-lumen portions of a dilatation catheter so as to maintain the inside diameters of the catheter lumens to as great an extent as possible.

Objects Of The Invention

It is one object of the present invention to provide an improved method of joining portions of an inflation/expansion catheter.

It is another object of the present invention to provide an improved method of joining portions of a dilatation catheter.

Summary Of The Invention

The objects of the present invention are accomplished by the method of joining portions of an inflation/expansion catheter as described in more detail below. In particular, the objects of the present invention will be described in relation to a dilatation catheter as shown in the drawing figures.

Brief Description Of The Drawings

- Fig. 1 is a plan view of a dilatation catheter according to a first embodiment of the present invention.
- Fig. 2 is a cross-sectional view of dilatation catheter according to the present invention, taken along line A-A of Fig. 1.
- Fig. 3 is a cross-sectional view of a dilatation catheter according to the present invention, taken along

3

line B-B of Fig. 1.

Fig. 4 is a plan view of the inner member of the balloon shaft of the dilatation catheter according to the present invention.

Fig. 5 is a plan view of the balloon member of the balloon shaft of the dilatation catheter according to the present invention.

Detailed Description Of The Invention

Fig. 1 is a plan view of a dilatation catheter according to a first embodiment of the present invention. The catheter, generally designated by reference numeral 10, includes a multi-lumen shaft 20, one lumen of which serves as a guidewire lumen, and the other lumen of which serves as a balloon inflation lumen. The multi-lumen shaft 20, is joined at connection portion 40, to two connection lumens 42, 46. The connection lumen 42, is connected to the guidewire lumen of the multi-lumen shaft 20, at the connection portion 40, and includes a hub connector 44, at the proximal end thereof. connection lumen 46, is connected to the inflation lumen of the multi-lumen shaft 20, at the connection portion 40, and includes a hub connector 48, at the proximal end thereof.

The guidewire lumen runs the entire length of the multi-lumen shaft 20, from a proximal end 24, to a joining portion 28. The inflation lumen also runs from the proximal end 24, of the multi-lumen shaft 20, to the joining portion 28. As shown most clearly in Fig. 2, the guidewire lumen 50, and inflation lumen 53, are in a

4

PCT/US94/04574

side-by-side configuration. As clearly shown, guidewire lumen 50, has a much larger cross-sectional area than that of the inflation lumen 53. lumen shaft 20, is connected to a balloon shaft generally designated by reference numeral 30, at the joining portion 28. As best shown in Fig. 3, balloon shaft 30, comprises two separate tubes arranged in a coaxial In particular, the balloon shaft 30, configuration. includes an inner member 35, having a guidewire lumen 34; and a balloon member 37, having an inflation lumen 36. The multi-lumen shaft 20, is connected to the balloon shaft 30, in such a manner that the guidewire lumen 50, of the multi-lumen shaft 20, and the guidewire lumen 34, of the inner member 35, form a continuous lumen running the entire length of the catheter 10; and the inflation lumen 53, of the multi-lumen shaft 20, and the inflation lumen 36, of balloon member 37, form a continuous lumen running from the proximal end of the catheter 10, to a balloon 38, formed along a portion of the balloon member 37. The method of the present invention for connecting the multi-lumen shaft 20, and the balloon shaft 30, is described in greater detail below.

The method of joining the multi-lumen shaft 20, and balloon shaft 30, according to the present invention comprises the following steps. Initially, the multi-lumen shaft 20, is cut to length. Tables 1 and 2 show appropriate lengths for the multi-lumen shaft 20, of varying dilatation balloon sizes. The cut on the distal end of the multi-lumen shaft 20, may be a diagonal cut to facilitate the attachment of the multi-lumen shaft 20, to the balloon shaft 30. Following cutting, the multi-lumen shaft 20, should be thoroughly cleaned on the outside with a lint-free wipe and freon, and on the inside by attachment to a compressed air source.

Fig. 4 is a plan view of the inner member 35, of the balloon shaft 30, according to the present invention. In particular, Fig. 4 shows placement of marker bands, preferably gold marker bands, 62, 64, on the inner member 35, of the balloon shaft 30. Tables 1 and 2 include lengths for the marker band section 60, of the inner member 35, i.e. measured form the outside of marker band 62, to the outside of marker band 64.

An inner member 35, prefitted with marker bands 62, 64, is cut to length according to the dimensions given in Tables 1 and 2. The inner member proximal shaft side 66, is measured from the outside of marker band 62, and the inner member distal shaft side 68, is measured from the outside of marker band 64. Following cutting, the inner member 35, should be thoroughly cleaned with a lint-free wipe and if necessary, degreased; for example, by using freen.

Two or three small mandrils, typical sizes of which are shown in Tables 1 and 2, are then inserted through the entire length of inflation lumen 53, of the multilumen shaft 20, such that a portion of the mandrils extend beyond each end of the multi-lumen shaft 20. The mandrils should extend 1 to 2 cm beyond the distal end of the multi-lumen shaft 20, to facilitate further processing. The sizes of the small mandrils should be chosen so as to ensure a snug fit between the mandrils and the inside diameter of the inflation lumen 53. Optionally, a pin gauge may be used to enlarge the guidewire lumen 50, of the multi-lumen shaft 20, to a distance of about 1 cm.

The inner member 35, of balloon shaft 30, is slidingly engaged over an appropriately sized mandril;

6

e.g. 0.021°, and is threaded over the mandril such that the mandril extends beyond both ends of the inner member 35. The inner member mandril is then inserted approximately 8 to 10 cm into the guidewire lumen 50, of the multi-lumen shaft 20. The inner member 35, is then threaded along the inner member mandril and the proximal shaft side 66, of the inner member 35, is inserted about 1 cm into the guidewire lumen 50, of the multi-lumen shaft 20.

The balloon member 37, of balloon shaft 30, has a length in accordance with the dimensions shown in Tables 3 and 4. As best shown in Fig. 5, the overall length of balloon member 37, is determined by addition; of the dimensions shown in Tables 3 and 4 for balloon 38, size; proximal shaft 72, length; and distal shaft 74, length. The balloon member 37, is slidingly engaged over the inner member mandril and the inner member 35, and is threaded over the mandril and the inner member 35, toward the distal end of the multi-lumen shaft 20. The balloon member 37, is then slidingly engaged over the two or three mandrils inserted into the inflation lumen 53, of multi-lumen shaft 20, and then over the distal end of the multi-lumen shaft 20, to create a lap joint of about 1 cm in length. The balloon member 37, may optionally include an adaptation expansion at the proximal end, so as to facilitate the engagement of the balloon member 37, over the multi-lumen shaft 20, and the creation of the lap joint.

A teflon tube, having the size indicated in Tables 3 and 4 is engaged over the multi-lumen shaft 20, and positioned to cover the lap joint. The teflon tube may include a short transverse slit at the distal end, to facilitate the positioning of the teflon tube over the

WO 94/25097 PCT/US94/04574

7

lap joint. A first glass capillary tube having the inside diameter shown in Tables 3 and 4 is then engaged over the multi-lumen shaft 20, and positioned to cover the teflon tube.

Various methods may be used to accomplish a good seal or bond between the multi-lumen shaft 20, and the balloon shaft 30, at the lap joint.

According to one method, when the teflon and glass capillary tubes are in place, the lap joint portion is placed in a heater block assembly set at a temperature and air flow indicated in Tables 3 and 4. The lap joint is positioned to be in the middle of the air flow nozzle and is rotated to achieve even heating. The lap joint is heated until it turns a white or milkish color. Care should be taken to heat both edges of the lap joint. Pressure should be applied to the lap joint during heating and subsequent cooling in air.

Following the heating and cooling steps, the glass capillary and teflon tubes are removed. The lap joint should be smooth and sealed. If not properly sealed, the heating and cooling steps should be repeated using the small second glass capillary tube indicated in Tables 3 and 4. When the lap joint is sealed properly, the mandrils may all be removed.

In a preferred method, the lap joint portion covered by the teflon and glass capillary tubes is placed in a heater coil and compression monitor assembly set at predetermined values for control of current and time. In particular, amps for the heater coil and compression force for the compression monitor may be set in accordance with the parameters indicated in Table 5. The

lap joint then cycles through a heating time, compression time, and cooling time, also as indicated in Table 5.

Following the heating, compression and cooling steps, the glass capillary and teflon tubes are removed. The lap joint should be smooth and sealed. If not properly sealed, the heating, compression and cooling steps may be repeated using the small second glass capillary tube indicated in Tables 3 and 4. When the lap joint is sealed properly, the mandrils may all be removed.

The provision of a diagonal cut on the multi-lumen shaft 20, provides several advantages in the above process. In particular, as noted above, the diagonal cut facilitates the attachment of the multi-lumen shaft 20, to the inner member 35, and the balloon member 37, of the balloon shaft 30. In addition, the diagonal cut reduces the amount of material which undergoes the sealing or bonding process. In particular, the amount of material located in the lap joint area is reduced, thereby making the bonding of the multi-lumen shaft 20, to the various components of the balloon shaft 30, easier to control and providing a more even and smooth bond area.

 \times \times \times \times \times

.025

.011

011

Ę

5 5 5 5

2.0 4.0 4.0

×

3.0

3.0

.045

×

.025

.011 .011 .011 .011

.011

5 5

109 109 109 109 109

109

4 6 4 6 1 6 4 4 4 4 4

8.6

थ थ थ थ थ थ थ

.045

.025

.045

.025

6

9

TABLE

Member Size Inner Lumen Distal Shaft Side (± 0.2 cm) Inner Member Inner Member Shaft Side (± 0.2 cm) Proximal Mandril 0.D. (inches) Marker Band mm) 0.5 Length Ħ Multi-Lumen E S 0 # Balloon Size CE CE × E E

and Shaft Length (Tables cm shorter than the Distal sure inner member is 1 cm ±0.5 Make adjustments if necessary. Make

.011

.011

.011

1.4

3.4

109

×

3.5

×

3.5

E 5 5

Make sure each end is rounded off 2.0~X 2.0 through 3.5 X 4.0; mandril length 125 cm ±5 cm. both mandrils before using.

2.5 X

2.5 X

2.0

2.0

TABLE 2

Inner Member Lumen Size	.025 X .045								
Inner Member Distal Shaft Side (± 0.2 cm)	Long ***								
Inner Member Proximal Shaft Side (± 0.2 cm)	*	* *	* *	* *	* *	* *	* *	* *	* *
Mandril 0.D. (inches)	.012 & .012	.012 & .012	.012 & .012	.012 & .012	.012 & .012	.012 & .012	.012 & .012	.012 & .012	.012 & .012
Marker Band Length (± 0.5 mm)	1.4 cm	3.4 cm	2.4 cm						
se Multi-Lumen	77 cm		77 cm	77 ст					
Balloon Size	_	4.0 X 4.0	5.0 X 2.0	5.0 X 4.0	6.0 x 2.0	6.0 x 4.0	7.0 X 2.0	7.0 X 4.0	8.0 x 3.0

Make sure each end is rounded off on 4.0 X 2.0 through 8.0 X 3.0; mandril length 90 cm ± 5 cm. both mandrils before using.

Cut inner member to fit balloon. Cut inner member even with Distal Shaft Length (Tables 3 and 4) after adaption is completed.

* *

SUBSTITUTE SHEET (RULE 26)

TABLE 3

Nozzle Air Flow (SCEH)	20-30	30-40	30-40	30-40	30-40	30-40	30-40	30-40
ure		250	270	270	270	270	270	270
stal Teflon 1st Glass 2nd Glass Temperat aft Size I.D. I.D. I.D. ingth (\pm 0.001°) (\pm .2 cm) (\pm .2 cm) (\pm 10 F)	.060/.061	.060/.061	.060/.061	.060/.061	.062/.063	.062/.063	.062/.063	.062/.063
1st Glass I.D. (± .2 cm)	.063/.064	.063/.064	.063/.064	.063/.064	.064/.065	.064/.065	.064/.065	.064/.065
Teflon Size	0.052	0.052	0.052	0.052	0.052	0.052	0.052	0.052
Distal Shaft Length (± 0.2 cm)	3 cm	3 cm	3 cm	3 cm	3 сш	3 ст	3 cm	3 cm
Proximal Shaft Length (± 2 mm)	7 cm							
Balloon Size (mm X cm)	2.0 X 2.0	2.0 X 4.0	2.5 x 2.0	2.5 X 4.0	3.0 x 2.0	3.0 X 4.0	3.5 x 2.0	3.5 X 4.0

Length of teflon is 5 cm ±5 mm.

Check temperature with probe to ensure proper temperature.

TABLE 4

Balloon Size	Proximal Shaft	Distal Shaft	Teflon Size	1st Glass I.D.	2nd Glass I.D.	Temperature	Nozzle Air Flow
(man X cm)	Length (± 2 mm)	Length (± 0.2 cm)	ength $(\pm 0.2 \text{ cm})$ $(\pm 0.001")$ $(\pm .2 \text{ cm})$ $(\pm .2 \text{ cm})$ $(\pm 10 \text{ F})$	(± .2 cm)	(± .2 cm)	(± 10 F)	(SCFH)
					i		
4.0 X 2.0	7 cm	2 cm	0.052	.066/.067	990./590.	280	30-40
4.0 X 4.0	7 cm	2 cm	0.052	.066/.067	990./590.	280	30-40
5.0 X 2.0	7 cm	2 cm	0.052	.067/.068	990./590.	280	30-40
5.0 X 4.0	7 cm	2 cm	0.052	.067/.068	990./590.	280	30-40
6.0 x 2.0	7 cm	2 cm	0.052	.067/.068	990./590.	290	30-40
6.0 x 4.0	7 cm	2 cm	0.052	.067/.068	990./590.	290	30-40
7.0 x 2.0	7 cm	2 cm	0.052	.067/.068	990./590.	290	30-40
7.0 X 4.0	7 cm .	2 cm	0.052	.067/.068	990./590.	290	30-40
8.0 X 3.0	7 cm	2 cm	0.052	690'/890'	.066/.067	290	30-40

Check temperature with probe to ensure proper temperature.

Length of teflon is 5 cm 15 mm.

SUBSTITUTE SHEET (RULE 26)

		TAB	TABLE 5		
Balloon Size (mm X cm)	Amps	Compression Force (Lbs)	Heating Time (Seconds)	Compression Time (Seconds)	Cooling Time (Seconds)
2.0 X 2.0	$6.30 \pm .15$.35 ± .15	26.0 ± .6	14.0 ± .6	10
2.0 X 4.0	$6.30 \pm .15$	$.35 \pm .15$	26.0 ± .6	14.0 ± .6	10
2.5 x 2.0	$6.30 \pm .15$.35 ± .15	$26.0 \pm .6$	14.0 ± .6	10
2.5 x 4.0	6.30 ± .15	.35 ± .15	26.0 ± .6	14.0 ± .6	10
3.0 x 2.0	$6.10 \pm .10$.40 ± .10	24.0 ± .4	14.0 ± .4	10
3.0 x 4.0	$6.10 \pm .10$.40 ± .10	24.0 ± .4	14.0 ± .4	10
3.5 x 2.0	$6.10 \pm .10$.40 ± .10	24.0 ± .4	14.0 ± .4	10
3.5 X 4.0	$6.10 \pm .10$.40 ± .10	24.0 ± .4	14.0 ± .4	10
4.0 x 2.0	$6.10 \pm .10$.40 ± .10	24.0 ± .4	14.0 ± .4	10
4.0 X 4.0	$6.10 \pm .10$.40 ± .10	24.0 ± .4	14.0 ± .4	10
5.0 x 2.0	6.10 ± .10	.40 ± .10	$24.0 \pm .4$	14.0 ± .4	10
5.0 X 4.0	$6.10 \pm .10$.40 ± .10	24.0 ± .4	14.0 ± .4	10
6.0 X 2.0	$6.10 \pm .05$.45 ± .15	26.0 ± .6	$16.0 \pm .6$	10
6.0 x 4.0	$6.15 \pm .10$.45 ± .15	26.0 ± .6	$16.0 \pm .6$	10
7.0 X 2.0	6.10 ± .10	.40 ± .10	$24.0 \pm .4$	14.0 ± .4	10
7.0 X 4.0	$6.10 \pm .10$.40 ± .10	24.0 ± .4	14.0 ± .4	10
8.0 x 3.0	$6.10 \pm .10$.40 ± .10	24.0 ± .4	14.0 ± .4	1.0

SUBSTITUTE SHEET (RULE 26)

By using the method according to the present invention, it is possible to securely and accurately join different portions of a dilatation catheter. In particular, the method of the present invention makes it possible to join a multi-lumen shaft to a coaxial shaft.

The method described above relating to a two lumen catheter may easily be extended to a catheter having any number of lumens.

Catheters having more than two lumens may be constructed in a very similar method to that described above for catheter 10. It will be easily recognized that the only differences will be the addition of further mandrils and the additional steps of connecting a third set of lumens at the joining portion 28.

In addition, while the above description relates specifically to a dilatation catheter, the method described could easily be used for joining separate portions of almost any type of catheter. For example, the catheter joined by the method according to the present invention may include a lumen for use in the delivery of drugs or other pharmacologic agents.

The foregoing has been a description of certain preferred embodiments of the present invention, but is not intended to limit the invention in any way. Rather, many modifications, variations and changes in details may be made within the scope of the present invention.

What is claimed is:

- A method of joining portions of a catheter, said catheter including
 - a proximal shaft with at least two lumens, wherein said at least two lumens are in a side-by-side arrangement;
 - a distal shaft with at least two tubes, wherein said at least two tubes are arranged in a coaxial arrangement;

said method comprising the steps of:

- a) connecting said tubes of said distal shaft to corresponding said lumens of said proximal shaft until only one tube remains for said distal shaft and only one lumen remains for said proximal shaft;
- b) connecting said remaining lumen of said proximal shaft to said remaining tube of said distal shaft to create a lap joint; and
- c) treating said lap joint to seal said proximal shaft to said distal shaft.
- 2. A method according to claim 1, wherein said proximal shaft includes a first lumen comprising a guidewire lumen, and a second lumen comprising an inflation lumen; and wherein said distal shaft includes a first tube comprising an inner member having a guidewire lumen, and a second tube comprising an outer member having an inflation lumen and having an inflation/expansion member formed along a portion of its length.

- 3. A method according to claim 2, wherein said step a) comprises connecting said first lumen of said proximal shaft to said first tube of said distal shaft, so as to form a continuous guidewire lumen; and wherein said step b) comprises connecting said second lumen of said proximal shaft and said second tube of said distal shaft, so as to form a continuous inflation lumen connected with said inflation/expansion member.
- 4. A method according to claim 2, wherein said catheter is a dilatation catheter and wherein said inflation/expansion member is a dilatation balloon.
- 5. A method according to claim 4, wherein said catheter is a percutaneous transluminal angioplasty catheter.
- 6. A method according to claim 1, wherein said proximal shaft includes a first lumen comprising a delivery lumen, and a second lumen comprising an inflation lumen; and wherein said distal shaft includes a first tube comprising an inner member having a delivery lumen, and a second tube comprising an outer member having an inflation lumen and having an inflation/expansion member formed along a portion of its length.

- 7. A method according to claim 6, wherein said step a) comprises connecting said first lumen of said proximal shaft to said first tube of said distal shaft, so as to form a continuous delivery lumen; and wherein said step b) comprises connecting said second lumen of said proximal shaft and said second tube of said distal shaft, so as to form a continuous inflation lumen connected with said inflation/expansion member.
- 8. A method according to claim 2, wherein said catheter is a drug delivery catheter and is capable of delivering drugs or other pharmacologic agents through said delivery lumen.
- 9. A method according to claim 1, wherein said proximal shaft includes only two lumens and said distal shaft includes only two tubes.
- 10. A method according to claim 1, wherein said proximal shaft includes more than two lumens and said distal shaft includes more than two tubes.

11. A method according to claim 1, wherein said step of treating said lap joint comprises:

positioning a teflon tube over the lap joint; positioning a glass capillary tube over the teflon tube to form an assembly;

placing said assembly within a heater block having an air flow nozzle, such that lap joint is positioned to be in the middle of the air flow nozzle;

rotating said assembly to achieve even heating of said lap joint;

cooling said assembly; and removing said glass capillary tube and said teflon tube.

12. A method according to claim 1, wherein said step of treating said lap joint comprises:

positioning a teflon tube over the lap joint; positioning a glass capillary tube over the teflon tube to form an assembly;

placing said assembly within a heater coil and compression monitor set at a predetermined current and time;

cycling said assembly through a heating time, compression time and cooling time; and

removing said glass capillary tube and said teflon tube.

13. A method of joining different types of tubes, said method comprising;

providing a first shaft with at least two lumens, wherein said at least two lumens are in a side-by-side arrangement;


providing a second shaft comprised of at least two tubes, wherein each tube has a lumen associated therewith, and wherein said at least two tubes are arranged in a coaxial arrangement;


connecting said tubes of said second shaft to corresponding said lumens of said first shaft until only one tube remains for said second shaft and only one lumen remains for said first shaft;


connecting said remaining lumen of said first shaft to said remaining tube of said second shaft to create a lap joint; and

treating said lap joint to seal said first shaft to said second shaft.

- 14. A method according to claim 13, wherein said first shaft includes only two lumens and second shaft includes only two tubes.
- 15. A method according to claim 13, wherein said first shaft includes more than two lumens and said second shaft includes more than two tubes.

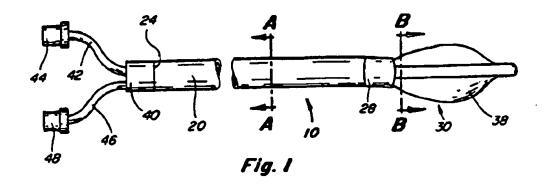

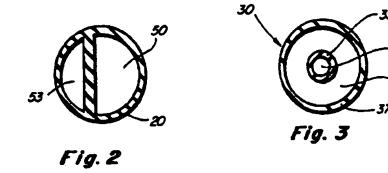

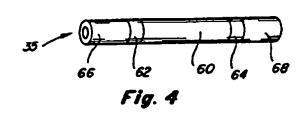

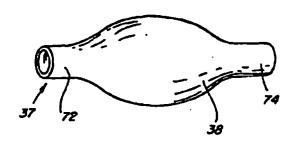


Fig. 5

INTERNATIONAL SEARCH REPORT

International application No. PCT/US94/04574

A. CLASSIFICATION OF SUBJECT MATTER							
IPC(5) :A61M 29/00 US CL :604/103							
According to International Patent Classification (IPC) or to both national classification and IPC							
B. FIELDS SEARCHED Minimum documentation searched (classification system followed by classification symbols)							
U.S. : 604/96-103; 606/192-194							
		to the California had					
Documentati NONE	ion searched other than minimum documentation to the extent that such documents are included	in the news serience					
NONE	ata base consulted during the international search (name of data base and, where practicable,	acaren erris esco)					
C. DOCUMENTS CONSIDERED TO BE RELEVANT							
Category*	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.					
X, P	US, A, 5,209,729, (HOFMANN ET AL.), 11 May 1993. See columns 2-4.	1-10, 13-15					
A, P	US, A, 5,295,961, (NIEDERHAUSER ET AL.), 22 March 1994.	1					
A, E	US, A, 5,306,247, (PFENNINGER), 26 April 1994.	1					
A	US, A, 5,002,559, (TOWER), 26 March 1991.See column 4, lines 4 and 5.	11, 12					
A	US, A, 4,913,701, (TOWER), 03 April 1990. See column 3, lines 48-68.	11, 12					
		·					
Further documents are listed in the communion of Box C. See patent family annex.							
Special categories of cited documents: The later document published after the international filing date or priority date and not in conflict with the application but sized to understand the							
to be of perticular reterence: the claimed invention cannot be							
E cartier decrement published on or after the immensional filling date considered sovel or cannot be considered to involve an inventive step when the document is taken alone.							
cized to antablish the publication date of another citation or other special reason (as specified) document of purticular relevance; the chained invention cannot be considered to involve an inventive seep when the document is							
O document referring to an oral disclosure, we, exhibition or other combined with one or more other such documents, such combination being obvious to a person skilled in the art							
·P· de	comment published prior to the interestional filing date but later than "&" document member of the same potent e priority date claimed						
	actual completion of the international search Date of mailing of the international search O 4 0	arch report CT 1994					
	mailing address of the ISA/US Authorized officer						
Commission Box PCT	oper of Patents and Trademarks FRANK WILKINS, III						
	nn, D.C. 20231 No. (703) 305-3230 Telephone No (703) 308-2681						

į